NovaLisa™

Mycoplasma pneumoniae

IgM – ELISA

Enzyme immunoassay for the qualitative determination of IgM-class antibodies against Mycoplasma pneumoniae in human serum or plasma

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For further languages please contact our authorized distributors.

Bibliography / Literatur / Bibliographie / Bibliografia / Bibliografía / Bibliografia	Page / Seite / Page / Pagina / Página/	34
Symbols Key / Symbolschlüssel / Explication des symboles / Legenda / Símbolos / Legenda dos Símbolos	Page / Seite / Page / Pagina / Página	35
Summary of Test Procedure/ Kurzanleitung Testdurchführung/ Résumé de la procédure de test/ Schema della procedura/ Resumen de la técnica / Resumo do Procedimento de Teste	Page / Seite / Page/ / Pagina / Página	36

Product Number: MYCM0350 (96 Determinations)

1. INTRODUCTION

The mycoplasmas belong to the class Mollicutes comprising three distinct families and four genera, one of which is *Mycoplasma* with over 60 species. Mycoplasmae are the smallest free-living organisms known (300 to 500 nm in diameter) and unlike regular bacteria they lack a cell wall. Mycoplasmas are extracellular parasites, especially on mucous membranes, which can cause infections in humans, animals, plants, and cell cultures. *Mycoplasma pneumoniae* is primarily a respiratory pathogen (obligate) in humans involving the nasopharynx, throat, trachea, bronchi, bronchioles, and alveoli. Other Mycoplasmae, *M. buccale*, *M. faucium*, *M. orale* and *M. salivarium* are commensals in the oral cavity. *Mycoplasma hominis* and *Ureaplasma urealyticum* inhabit primarily the genital tract and may act as opportunistic invaders. *M. pneumoniae* is by far the most important pathogen of this group. Infection with *M. pneumoniae* occurs worldwide, its epidemiology has been studied primarily in the USA, Europe, and Japan. Infections are endemic in larger urban areas, and epidemic increases are observed at varying intervals. *Mycoplasma pneumoniae* has been estimated to cause 15-20% of all pneumonias; the rate is highest in children and young adults. 74% of infections with *M. pneumoniae* are asymptomatic, reinfection may occur. Naturally acquired immunity to infection with *M. pneumoniae* appears to be of limited duration (2-3 years).

Species	Disease	Mechanism of Infection
<i>Mycoplasma pneumoniae</i>	Respiratory tract disease: From asymptomatic infection to pharyngitis, bronchitis, croup, tracheobronchitis, pneumonitis and pneumonia	It is not known whether spread is primarily by droplets, direct or indirect contact, or by all these means

Infection may be identified by:

- Microscopy
- Hemadsorption, Tetrazolium reduction (specific for *M. pneumoniae*), Complement fixation
- Detection of antibody production by ELISA: Serology: complement fixation (CF), neutralization (N) and hemagglutination-inhibition (HAI); Detection of antibodies and the hexon antigen by ELISA.

2. INTENDED USE

The NovaTec *Mycoplasma pneumoniae* IgM-ELISA is intended for the qualitative determination of IgM class antibodies against *M. pneumoniae* in human serum or plasma (citrate).

3. PRINCIPLE OF THE ASSAY

The qualitative immunoenzymatic determination of IgM-class antibodies against *Mycoplasma pneumoniae* is based on the ELISA (Enzyme-linked Immunosorbent Assay) technique.

Microtiter strip wells are precoated with *Mycoplasma pneumoniae* antigens to bind corresponding antibodies of the specimen. After washing the wells to remove all unbound sample material horseradish peroxidase (HRP) labelled anti-human IgM conjugate is added. This conjugate binds to the captured *Mycoplasma pneumoniae*-specific antibodies. The immune complex formed by the bound conjugate is visualized by adding Tetramethylbenzidine (TMB) substrate which gives a blue reaction product. The intensity of this product is proportional to the amount of *Mycoplasma pneumoniae*-specific IgM antibodies in the specimen. Sulphuric acid is added to stop the reaction. This produces a yellow endpoint colour. Absorbance at 450 nm is read using an ELISA microwell plate reader.

4. MATERIALS

4.1. Reagents supplied

- ***Mycoplasma pneumoniae* Coated Wells (IgM):** 12 breakapart 8-well snap-off strips coated with *Mycoplasma pneumoniae* antigen; in resealable aluminium foil.
- **IgM Sample Diluent ***:** 1 bottle containing 100 ml of buffer for sample dilution; pH 7.2 ± 0.2; coloured green; ready to use; white cap.
- **Stop Solution:** 1 bottle containing 15 ml sulphuric acid, 0.2 mol/l; ready to use; red cap.
- **Washing Solution (20x conc.):*** 1 bottle containing 50 ml of a 20-fold concentrated buffer (pH 7.2 ± 0.2) for washing the wells; white cap.
- ***Mycoplasma pneumoniae* anti-IgM Conjugate**:** 1 bottle containing 20 ml of peroxidase labelled rabbit antibody to human IgM; coloured red, ready to use; black cap.
- **TMB Substrate Solution:** 1 bottle containing 15 ml 3,3',5,5'-tetramethylbenzidine (TMB); ready to use; yellow cap.
- ***Mycoplasma pneumoniae* IgM Positive Control***:** 1 bottle containing 2 ml; coloured yellow; ready to use; red cap.
- ***Mycoplasma pneumoniae* IgM Cut-off Control***:** 1 bottle containing 3 ml; coloured yellow; ready to use; green cap.
- ***Mycoplasma pneumoniae* IgM Negative Control***:** 1 bottle containing 2 ml; coloured yellow; ready to use; blue cap.

* contains 0.1 % Bronidox L after dilution

** contains 0.2 % Bronidox L

*** contains 0.1 % Kathon

4.2. Materials supplied

- 1 Strip holder
- 1 Cover foil
- 1 Test protocol
- 1 distribution and identification plan

4.3. Materials and Equipment needed

- ELISA microwell plate reader, equipped for the measurement of absorbance at 450/620 nm
- Incubator 37°C
- Manual or automatic equipment for rinsing wells
- Pipettes to deliver volumes between 10 and 1000 µl
- Vortex tube mixer
- Deionised or (freshly) distilled water
- Disposable tubes
- Timer

5. STABILITY AND STORAGE

The reagents are stable up to the expiry date stated on the label when stored at 2...8 °C.

6. REAGENT PREPARATION

It is very important to bring all reagents, samples and controls to room temperature (20...25°C) before starting the test run!

6.1. Coated snap-off Strips

The ready to use breakapart snap-off strips are coated with Mycoplasma pneumoniae antigen. Store at 2...8°C. *Immediately after removal of strips, the remaining strips should be resealed in the aluminium foil along with the desiccant supplied and stored at 2...8 °C; stability until expiry date.*

6.2. Mycoplasma pneumoniae anti-IgM Conjugate

The bottle contains 20 ml of a solution with anti-human-IgM horseradish peroxidase, buffer, stabilizers, preservatives and an inert red dye. The solution is ready to use. Store at 2...8°C. *After first opening stability until expiry date when stored at 2...8°C.*

6.3. Controls

The bottles labelled with Positive, Cut-off and Negative Control contain a ready to use control solution. It contains 0.1% Kathon and has to be stored at 2...8°C. *After first opening stability until expiry date when stored at 2...8°C.*

6.4. IgM Sample Diluent

The bottle contains 100 ml phosphate buffer, anti human-IgG, stabilizers, preservatives and an inert green dye. It is used for the dilution of the patient specimen. The solution contains antihuman IgG class antibodies to eliminate competitive inhibition from specific IgG class antibody to remove rheumatoid factor. This ready to use solution has to be stored at 2...8°C. *After first opening stability until expiry date when stored at 2...8°C.*

6.5. Washing Solution (20xconc.)

The bottle contains 50 ml of a concentrated buffer, detergents and preservatives. Dilute washing solution 1+19; e.g. 10 ml washing solution + 190 ml fresh and germ free redistilled water. The diluted buffer is stable for 5 days at room temperature. *Crystals in the solution disappear by warming up to 37 °C in a water bath. After first opening the concentrate is stable until the expiry date.*

6.6. TMB Substrate Solution

The bottle contains 15 ml of a tetramethylbenzidine/hydrogen peroxide system. The reagent is ready to use and has to be stored at 2...8°C, away from the light. *The solution should be colourless or could have a slight blue tinge. If the substrate turns into blue, it may have become contaminated and should be thrown away. After first opening stability until expiry date when stored at 2...8°C.*

6.7. Stop Solution

The bottle contains 15 ml 0.2 M sulphuric acid solution (R 36/38, S 26). This ready to use solution has to be stored at 2...8°C. *After first opening stability until expiry date.*

7. SPECIMEN COLLECTION AND PREPARATION

Use human serum or plasma (citrate) samples with this assay. If the assay is performed within 5 days after sample collection, the specimen should be kept at 2...8°C; otherwise they should be aliquoted and stored deep-frozen (-20 to -70°C). If samples are stored frozen, mix thawed samples well before testing. *Avoid repeated freezing and thawing.*
Heat inactivation of samples is not recommended.

7.1. Sample Dilution

Before assaying, all samples should be diluted 1+100 with IgM Sample Diluent. Dispense 10µl sample and 1 ml IgM Sample Diluent into tubes to obtain a 1+100 dilution and thoroughly mix with a Vortex.

8. ASSAY PROCEDURE

8.1. Test Preparation

Please read the test protocol carefully **before** performing the assay. Result reliability depends on strict adherence to the test protocol as described. The following test procedure is only validated for manual procedure. If performing the test on ELISA automatic systems we recommend to increase the washing steps from three to five and the volume of washing solution from 300µl to 350µl to avoid washing effects. Prior to commencing the assay, the distribution and identification plan for all specimens and controls should be carefully established on the result sheet supplied in the kit. Select the required number of microtiter strips or wells and insert them into the holder.

Please allocate at least:

1 well	(e.g. A1)	for the substrate blank,
1 well	(e.g. B1)	for the negative control,
2 wells	(e.g. C1+D1)	for the cut-off control and
1 well	(e.g. E1)	for the positive control.

It is recommended to determine controls and patient samples in duplicate, if necessary.

Perform all assay steps in the order given and without any appreciable delays between the steps.

A clean, disposable tip should be used for dispensing each control and sample.

Adjust the incubator to $37^{\circ} \pm 1^{\circ}\text{C}$.

1. Dispense 100µl controls and diluted samples into their respective wells. Leave well A1 for substrate blank.
2. Cover wells with the foil supplied in the kit.
3. **Incubate for 1 hour \pm 5 min at $37 \pm 1^{\circ}\text{C}$.**
4. When incubation has been completed, remove the foil, aspirate the content of the wells and wash each well three times with 300µl of Washing Solution. Avoid overflows from the reaction wells. The soak time between each wash cycle should be >5 sec. At the end carefully remove remaining fluid by tapping strips on tissue paper prior to the next step!
Note: Washing is critical! Insufficient washing results in poor precision and falsely elevated absorbance values.
5. Dispense 100µl Mycoplasma pneumoniae anti-IgM Conjugate into all wells except for the blank well (e.g. A1). Cover with foil.
6. **Incubate for 30 min at room temperature. Do not expose to direct sunlight.**
7. Repeat step 4.
8. Dispense 100µl TMB Substrate Solution into all wells
9. **Incubate for exactly 15 min at room temperature in the dark.**
10. Dispense 100µl Stop Solution into all wells in the same order and at the same rate as for the TMB Substrate Solution.
Any blue colour developed during the incubation turns into yellow.
Note: Highly positive patient samples can cause dark precipitates of the chromogen! These precipitates have an influence when reading the optical density. Predilution of the sample with physiological sodium chloride solution, for example 1+1, is recommended. Then dilute the sample 1+100 with dilution buffer and multiply the results in NTU by 2.
11. Measure the absorbance of the specimen at 450/620 nm within 30 min after addition of the Stop Solution.

8.2. Measurement

Adjust the ELISA Microwell Plate Reader **to zero** using the **substrate blank in well A1**.

If - due to technical reasons - the ELISA reader cannot be adjusted to zero using the substrate blank in well A1, subtract the absorbance value of well A1 from all other absorbance values measured in order to obtain reliable results!

Measure the absorbance of all wells at **450 nm** and record the absorbance values for each control and patient sample in the distribution and identification plan.

Dual wavelength reading using 620 nm as reference wavelength is recommended.

Where applicable calculate the **mean absorbance values** of all duplicates.

9. RESULTS

9.1. Run Validation Criteria

In order for an assay to be considered valid, the following criteria must be met:

- **Substrate blank** in A1: Absorbance value < **0.100**.
- **Negative control** in B1: Absorbance value < **0.200 and < cut-off**
- **Cut-off control** in C1 and D1: Absorbance value **0.150 – 1.30**.
- **Positive control** in E1: Absorbance value > **cut-off**.

If these criteria are not met, the test is not valid and must be repeated.

9.2. Calculation of Results

The cut-off is the mean absorbance value of the Cut-off control determinations.

Example: Absorbance value Cut-off control 0.49 + absorbance value Cut-off control 0.47 = 0.96 / 2 = 0.48

$$\text{Cut-off} = 0.48$$

9.3. Interpretation of Results

Samples are considered **POSITIVE** if the absorbance value is higher than 10% over the cut-off.

Samples with an absorbance value of 10% above or below the cut-off should not be considered as clearly positive or negative
→ **grey zone**

It is recommended to repeat the test again 2 - 4 weeks later with a fresh sample. If results in the second test are again in the grey zone the sample has to be considered **NEGATIVE**.

Samples are considered **NEGATIVE** if the absorbance value is lower than 10% below the cut-off.

9.3.1. Results in NovaTec Units

Patient (mean) absorbance value x 10 / Cut-off = [NovaTec-Units = NTU]

Example: $\frac{1.786 \times 10}{0.48} = 37 \text{ NTU (NovaTec Units)}$

Cut-off:	10	NTU
Grey zone:	9-11	NTU
Negative:	<9	NTU
Positive:	>11	NTU

10. SPECIFIC PERFORMANCE CHARACTERISTICS

10.1. Precision

Interassay	n	Mean (NTU)	Cv (%)
Pos. serum	12	21.6	4.8
Pos. serum	12	13.3	7.8
Intraassay	n	Mean (E)	Cv (%)
Pos. serum	31	0.99	2.9
Pos. serum	22	0.76	9.2

10.2. Diagnostic specificity

The diagnostic specificity is defined as the probability of the assay of scoring negative in the absence of the specific analyte. It is >95 %.

10.3. Diagnostic sensitivity

The diagnostic sensitivity is defined as the probability of the assay of scoring positive in the presence of the specific analyte. It is 94.4 %.

10.4. Interferences

Interferences with hemolytic, lipemic or icteric sera are not observed up to a concentration of 10 mg/ml hemoglobin, 5 mg/ml triglycerides and 0.2 mg/ml bilirubin.

Note: The results refer to the groups of samples investigated; these are not guaranteed specifications.
--

11. LIMITATIONS OF THE PROCEDURE

Bacterial contamination or repeated freeze-thaw cycles of the specimen may affect the absorbance values. Diagnosis of an infectious disease should not be established on the basis of a single test result. A precise diagnosis should take into consideration clinical history, symptomatology as well as serological data.

In immunocompromised patients and newborns serological data only have restricted value.

12. PRECAUTIONS AND WARNINGS

- In compliance with article 1 paragraph 2b European directive 98/79/EC the use of the in vitro diagnostic medical devices is intended by the manufacturer to secure suitability, performances and safety of the product. Therefore the test procedure, the information, the precautions and warnings in the instructions for use have to be strictly followed. The use of the testkits with analyzers and similar equipment has to be validated. Any change in design, composition and test procedure as well as for any use in combination with other products not approved by the manufacturer is not authorized; the user himself is responsible for such changes. The manufacturer is not liable for false results and incidents for these reasons. The manufacturer is not liable for any results by visual analysis of the patient samples.
- Only for in-vitro diagnostic use.
- All components of human origin used for the production of these reagents have been tested for anti-HIV antibodies, anti-HCV antibodies and HBsAg and have been found to be non-reactive. Nevertheless, all materials should still be regarded and handled as potentially infectious.
- Do not interchange reagents or strips of different production lots.
- No reagents of other manufacturers should be used along with reagents of this test kit.
- Do not use reagents after expiry date stated on the label.
- Use only clean pipette tips, dispensers, and lab ware.
- Do not interchange screw caps of reagent vials to avoid cross-contamination.
- Close reagent vials tightly immediately after use to avoid evaporation and microbial contamination.
- After first opening and subsequent storage check conjugate and control vials for microbial contamination prior to further use.
- To avoid cross-contamination and falsely elevated results pipette patient samples and dispense conjugate without splashing accurately to the bottom of wells.
- The NovaLisa™ ELISA is only designed for qualified personnel who are familiar with good laboratory practice.

WARNING:	In the used concentration Bronidox L has hardly any toxicological risk upon contact with skin and mucous membranes!
----------	---

WARNING:	Sulphuric acid irritates eyes and skin. Keep out of the reach of children. Upon contact with the eyes, rinse thoroughly with water and consult a doctor!
----------	--

12.1. Disposal Considerations

Residues of chemicals and preparations are generally considered as hazardous waste. The disposal of this kind of waste is regulated through national and regional laws and regulations. Contact your local authorities or waste management companies which will give advice on how to dispose hazardous waste.

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Clyde WA Jr: Clinical overview of typical *M. pneumoniae* infections. *Clin. Infect. Dis* 1993 Aug; 17 Suppl 1: S32-6 Medline

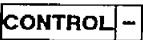
E. Jacobs. 1993. Serological Diagnosis of *M. pneumoniae* Infections: Current Procedures. *Clinical Infectious Diseases*. 17(suppl 1) 79-82

Foy, H.M. 1993. Infections caused by *M. pneumoniae* and possible populations of patients. *Clinical Infectious Diseases*. 17 (suppl 1) 37-46

Handley J.G., and L.D. Gray. 1997. The incidence of *M. pneumoniae* in Board Family Practice 11(6): 425-429

Vikerfors, T., Brodin, G., Grandien, M., Hirshberg, L., Krook, A., and C.A. Pettersson. Specific IgM antibodies for the diagnosis of *M. pneumoniae* infections. *Scand. J. Infect. Dis.* 20(6): 601-610

Wiegand, R.: *M. pneumoniae*. In: Brandis, H., W. Köhler, H.J. Eggers, G. Pulverer (ed.): *Med. Mikrobiologie* 1994. Gustav Fischer Verlag Stuttgart, Jena, New York (1994) 767-769

Symbols Key/ Symbolschlüssel/ Explication des symboles / Legenda / Símbolos / Legenda dos Símbolos	
	Manufactured by / Hergestellt von/ Fabriqué par/ Prodotto da/ Fabricado por / Fabricado por
	In Vitro Diagnostic Medical Device/ In Vitro Diagnosticum/ Dispositif médical de diagnostic <i>in vitro</i> / Diganostico <i>in vitro</i> / Producto para diagnóstico In vitro / Dispositivo Médico para Diagnóstico In Vitro
	Lot Number/ Chargenbezeichnung/ Numéro de lot/ Lotto/ Número de lote / Número de lote
	Expiration Date/ Verfallsdatum/ Date de péremption/ Scadenza/ Fecha de caducidad / Data de Validade
	Storage Temperature/ Lagertemperatur/ Température de conservation/ Temperatura di conservazione / Temperatura de almacenamiento / Temperatura de Armazenamento
	CE Mark/ CE-Zeichen/ Marquage CE / Marchio CE/ Marca CE / Marca CE
	Catalogue Number/ Katalog Nummer/ Référence du catalogue/ Numero di codice/ Número de Catálogo / Número de Catálogo
	Consult Instructions for Use/ Gebrauchsanweisung beachten/ Consulter la notice d'utilisation/ Consultare le istruzioni/ Consulte las Instrucciones de Uso / Consultar as Instruções de Utilização
MTP	Microplate/ Mikrotiterplatte/ Microplaque/ Micropiastra/ Microplaca. / Microplaca
CONJ	Conjugate/ Konjugat/ Conjugué/ Coniugato/ Conjugado / Conjugado
	Control serum, negative/ Kontrollserum, negative/ Sérum de contrôle négatif/ siero di controllo, negativo / Suero control negativo/ Soro de controle negativo / Soro de controllo negativo
	Control serum, positive/ Kontrollserum, positiv/ Sérum de contrôle positif/ siero di controllo, positivo/ Suero de control positivo / Soro de controle positivo / Soro de controllo positivo
	Cut off control serum/ Cut off Kontrollserum/ Sérum de contrôle du cut-off/ siero di controllo, cut-off/ Suero control Cut-off / Soro de controle Cut-off
DIL G	Sample diluent buffer IgM/ IgM-Probenverdünnungspuffer/ Tampon diluant pour échantillon IgM/ soluzione tampone per i campioni IgM/ solución tampón para muestras IgM / Tampão diluente para amostras IgM
SOLN STOP	Stop solution/ Stopplösung/ Solution d'arrêt/Soluzione bloccante / Solução de paragem
SUB TMB	TMB Substrate solution/ TMB-Substratlösung/ Substrat TMB/ soluzione substrato TMB/ solución substrato TMB / Solução substrato TMB
WASHBUF 20x	Washing solution 20x concentrated/ Waschlösung 20x konzentriert/ Solution de lavage concentré 20 x/ soluzione di lavaggio concentrazione x20/ solución de lavado concentrado x20 / Solução de lavagem concentrada 20x
	Contains sufficient for "n" tests/ Ausreichend für "n" Tests/ Contenu suffisant pour "n" tests/ Contenuto sufficiente per "n" saggi/ Contenido suficiente para "n" tests / Conteúdo suficiente para "n" testes

SCHEME OF THE ASSAY

Mycoplasma pneumoniae IgM-ELISA

Test Preparation

Prepare reagents and samples as described.
Establish the distribution and identification plan for all specimens and controls on the result sheet supplied in the kit.
Select the required number of microtiter strips or wells and insert them into the holder.

Assay Procedure

	Substrate blank (e.g. A1)	Negative control	Positive control	Cut-off control	Sample (diluted 1+100)
Negative control	-	100µl	-	-	-
Positive control	-	-	100µl	-	-
Cut-off control	-	-	-	100µl	-
Sample (diluted 1+100)	-	-	-	-	100µl
Cover wells with foil supplied in the kit Incubate for 1 h at 37°C Wash each well three times with 300µl of washing solution					
Conjugate	-	100µl	100µl	100µl	100µl
Cover wells with foil supplied in the kit Incubate for 30 min at room temperature Wash each well three times with 300µl of washing solution					
TMB Substrate	100µl	100µl	100µl	100µl	100µl
Incubate for exactly 15 min at room temperature in the dark					
Stop Solution	100µl	100µl	100µl	100µl	100µl
Photometric measurement at 450 nm (reference wavelength: 620 nm)					

NovaTec Immundiagnostica GmbH

Technologie & Waldpark

Waldstr. 23 A6
D-63128 Dietzenbach, Germany

Tel.: +49 (0) 6074-48760 Fax: +49 (0) 6074-487629

Email : info@NovaTec-ID.com

Internet: www.NovaTec-ID.com

MYCM0350engl,dt,fr,it,es,port-23072010-CR